

### Dra III



Product No : RE1256  
 Quantity : 200u



Lot :  
 Expiry Date :  
 Concentration : 5u/μl  
 Supplied with : 1ml of 10X Buffer *Dra III*  
 1ml of 10X Buffer UB  
 0.5ml Diluent Viva Buffer A  
 (BSA included in all Reaction Buffer)

Store at -20°C



info@vivantechnologies.com

#### Reaction Conditions:

**Buffer *Dra III*** ,  
 10mM Tris-HCl (pH 7.6), 10mM MgCl<sub>2</sub>,  
 200mM KCl, and 100μg/ml BSA.

**Incubate at 37°C.**

#### Dilution: Viva Buffer A

10mM Tris-HCl (pH 7.4 at 25°C), 50mM KCl, 0.1mM EDTA,  
 1mM DTT, 200μg/ml BSA and 50% glycerol.

**Thermal Inactivation:** 65°C for 20 minutes

#### Storage Buffer:

10mM Tris-HCl (pH 7.5), 200mM KCl, 0.1mM EDTA,  
 7mM 2-mercaptoethanol, 100μg/ml BSA and 50% glycerol.

#### Unit Definition:

1u is defined as the amount of enzyme that is required to digest 1μg of DNA in 1 hour at 37°C in 50μl of assay buffer.

#### Quality Control Assays:


##### Ligation/ Recutting Assay:

After 5-fold overdigestion with ***Dra III***, 70% of the DNA fragments can be ligated and recut.

##### Overdigestion assay:

An unaltered banding pattern was observed after 1μg of DNA was digested with 10u of ***Dra III*** for 16 hours at 37°C.

λ DNA  
 1.0% Agarose

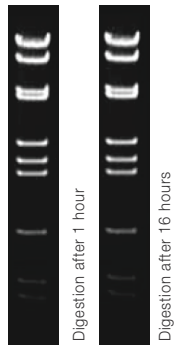

Activity in Reaction Buffer				
V1	V2	V3	V4	V5
75%	50%	50%	75%	75%

Buffer UB			
0.5X	1.0X	1.5X	2.0X
75%	75%	50%	50%

\* Buffer UB is provided for double digestion purpose.

#### NOTE:

- \* High enzyme concentration may result in **Star Activity**.
- \* Total reaction volume dependent on experiment.
- \* The amount of enzyme to be used is very much dependent on the DNA template.
- \* For plasmid DNA, 5-10X more enzyme is required.

#### Example of Digestion Reaction

Enzyme : 1 unit  
 Lambda 0.3μg/μl : 3.33μl (1μg DNA)  
 5X Reaction Buffer : 10μl  
 Sterile Distilled Water : Up to 50μl

Product Use Limitation

This product is for research purposes and *in vitro* use only.